




DNA離心濃縮機

DNA CONCENTRATOR CES-8000 + VP-60D

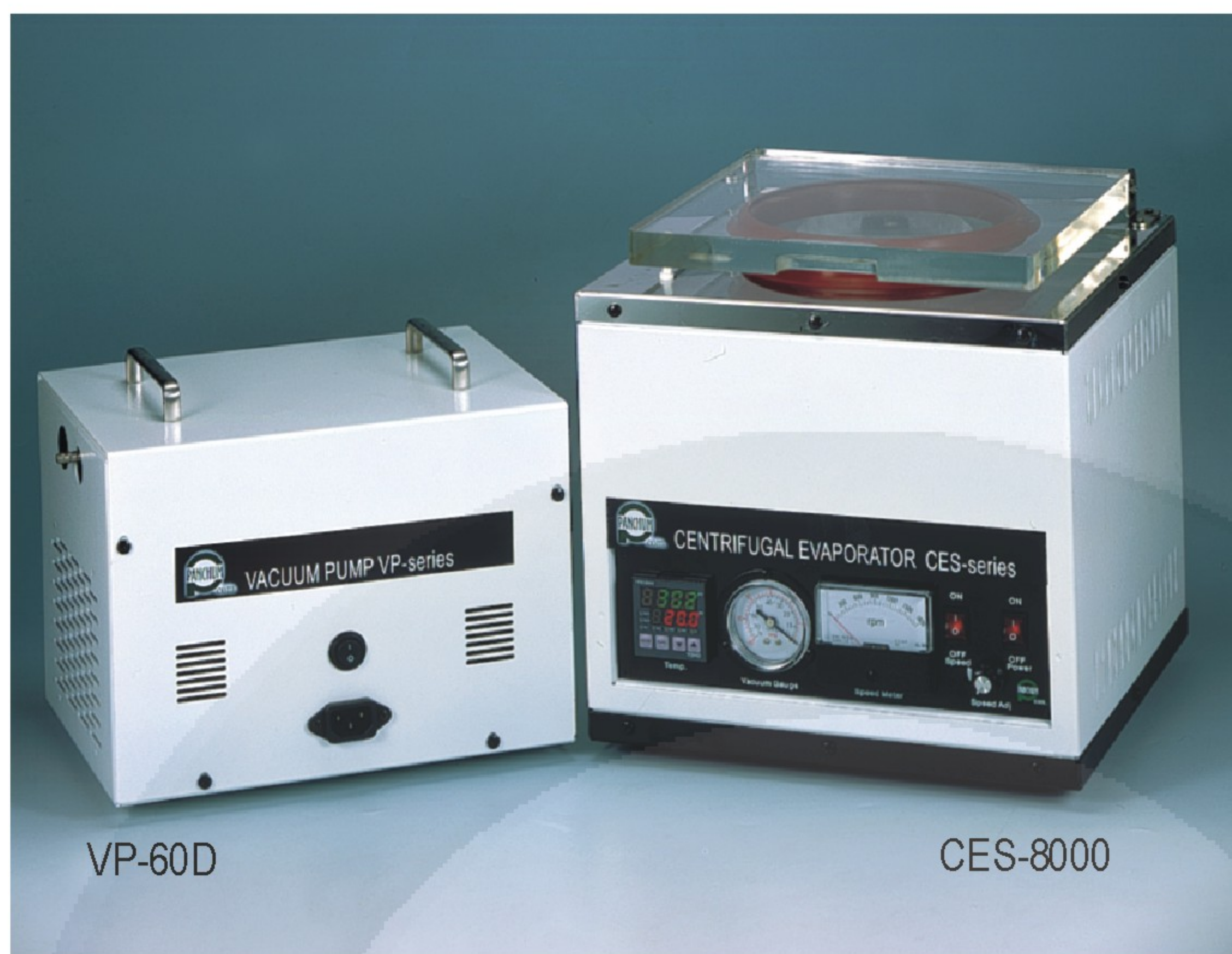

特點:

- 結合真空, 離心, 加熱於一體.
可避免樣品的突沸, 發泡, 並可快速濃縮.
- 方便監控溫度及真空度.
- 磁偶傳動式, 可獲得最佳真空效果.

SPECIFICATIONS

| | |
|-------------------------|-------------------------------------|
| MODEL | DNA CONCENTRATOR CES-8000+VP-60D |
| ROTOR SPEED | 100~1,800 rpm with meter |
| TEMP. RANGE | R.T.+5°C ~80°C |
| TEMP. STABILITY | ±0.1°C |
| TEMP. SETTING & DISPLAY | LED display |
| DRIVE SYSTEM | Magnetically coupled |
| VACUUM GAUGE | 0~760mmHg |
| CHAMBER MATERIAL | SUS304 plus teflon coating |
| CHAMBER DIMENSIONS | φ 142 x 110mm H |
| ROTOR | Optional |
| TOTAL DIMENSIONS(WxDxH) | 330 x 300 x 290mm |
| POWER SOURCE | 110/220V, 50/60Hz |
| VACUUM PUMP | Oilless type, 60l/min., 10 Torr |
| PRESSURE-SIDE TRAP | 500c.c. Bottle, optional |

** (本公司保有修改規格之權利, 不另行通知。)

We reserve the rights to change the specifications without notice.)

Optional rotor:

R-8000A, 1.5ml x 12 pcs

R-8000B, 1.5ml x 24 pcs

R-8000C, φ13 x 100mmL x 8 pcs or φ12 x 57mmL x 16 pcs

Specified for scintillation vial is available